

Transformation of MtrCAB and CysDes

Recorder: Yitian Zhou

Date: Jun. 13 2018

Introduction

Transformation of plasmid pSB1C3 containing MtrACB and CysDes to top10 *E.coil* strain.


由用户输入的一段描述性文字，介绍实验内容/目的。

Materials and equipment

由 Experiment Procedure 中每个小实验的标签自动生成的列表，包含实验所需耗材和仪器

Outline

由 basic experiment 自动生成的流程图，包含时间

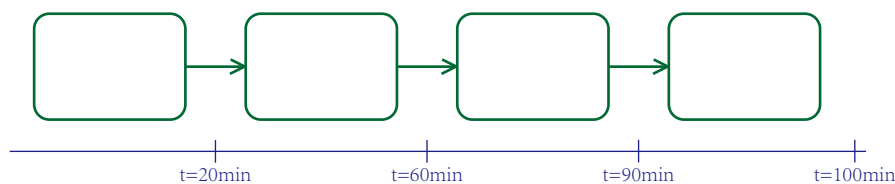

Result

由用户输入的实验结果，需要添加图片/表格/附件功能

Experiment Procedure

- 1 Take the competent bacteria from refrigerator and incubate them into ice about 5 mins until it is dissolved
- 2 Absorb 1 to 2 uL plasmid and mix it with bacteria solution thoroughly.
- 3 Put the tubes on the ice about 30 mins.
- 4 Make a heat shock at 42 degree centigrade about 60 sec
- 5 Put the tubes on the ice about 2 to 3 mins again.
- 6 Add 900 ul LB medium into EP tubes and cultivate the bacteria at 37 degree centigrade about 40 to 60 min.
- 7 Centrifuge them at 4,000 rpm about 2 min and we will see sediment in the tubes.
- 8 Discard the supernatant liquid and leave 220 ul medium.
- 9 Coat plate: add 200 ul solution in a large plate while add 20 ul solution in a small plate.
- 10 Cultivate these bacteria overnight for further use.

Notes

DO NOT add more than 5% volumn of bacteria solution
operate this step tenderly!

(Time SHOULD BE ACCURATE)

(Time SHOULD BE ACCURATE)

TIME COST about 1h

But we only have one flat plate (PET22 in AP plate) grows better. We can't get single colony in other plates. May be the heat-shock time we used isn't probable.

使用不同颜色标注

如：用户输入的关键参数用红色标出，用户添加的注释用绿色，系统注释用蓝色